


## LIFE SCIENCES

Nanoscribe's Photonic Professional GT 3D printers allow to easily produce tailored and reproducible 3D scaffolds and matrices. Here, the tailoring is achieved by on one hand a wide range of processable materials and on the other hand a resolution and feature sizes which, in combination, allow to mimic in vivo physiological environments. Accessible patterning scales cover sub-micrometer features and sample volumes of up to several 10 mm<sup>3</sup> at the same time. Samples can be produced from biocompatible, cell binding/repelling, hydrophobic and hydrophilic polymers with different mechanical characteristics as well as hydrogels.

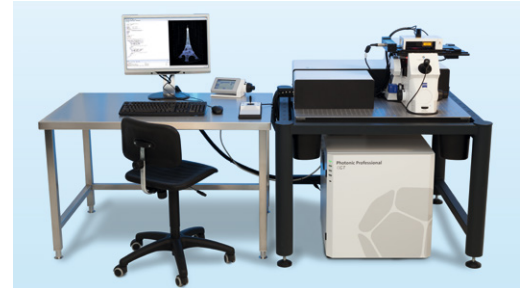


### CELL / EXTRACELLULAR MATRIX INTERACTION

**Challenge:** Fabrication of a microstructured 3D cell scaffold consisting of a cell-repellent material (grey) and a cell adhesion promoting material (red) as per design.

**Solution:** First, a core scaffold made of PEG-DA mixed with PETA (grey) is additively printed. Then, a second printing step adds cell adhesion points made of OrmoComp (red) precisely at the designed locations.

**Source:** DOI: 10.1002/adma.201004060

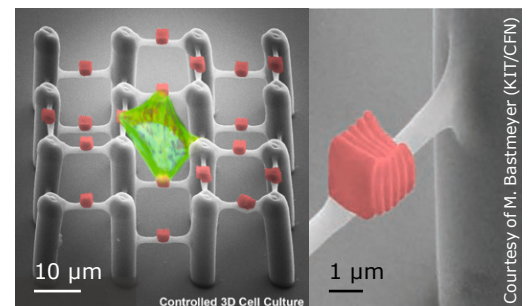


### TARGETED DRUG DELIVERY

**Challenge:** One-step fabrication of cargo-loaded helical swimming microrobots made of biodegradable superparamagnetic hydrogel for minimal invasive intervention.

**Solution:** A hydrogel consisting of PEG-DA/PETA, magnetite nanoparticles and photoinitiator is 3D printed. For loading, the helices are immersed in methylene blue.

**Source:** DOI: 10.1002/adma.201503112

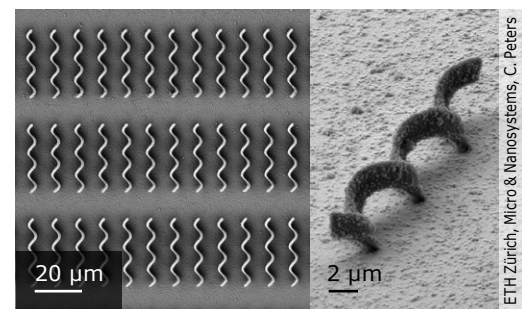


### CELL MIGRATION

**Challenge:** Fabrication of cage-like 3D micro-structures with sufficiently small pore size to differentiate tumorigenic from non-tumorigenic human breast cancer cells based on different invasion behavior.

**Solution:** Cubic cages with meshes of approximately 3, 18 and 85 µm<sup>2</sup> pore size as side faces are printed in Nanoscribe's biocompatible photoresist IP-L 780.

**Source:** DOI: 10.1038/srep10531

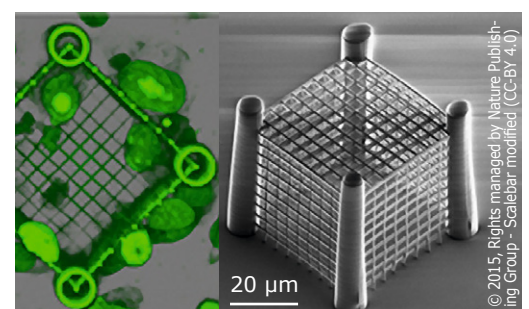


### TISSUE ENGINEERING / 3D CELL CULTURE

**Challenge:** High-fidelity replication of a human trabecular bone structure from high resolution 3D data obtained from a µ-CT scan.

**Solution:** µ-CT data is transferred to the standard 3D printing workflow. The scanned 3D object is reproduced with sub-micrometer precision in OrmoComp.

**Source:** DOI: 10.1016/j.actbio.2014.05.032

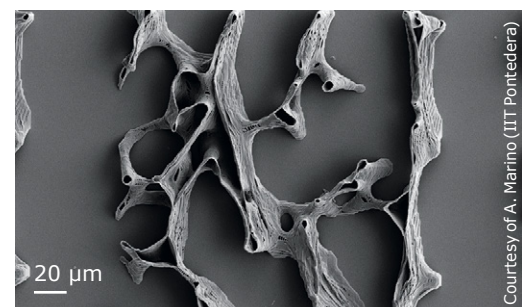

Learn more about the world's  
highest resolution 3D printer  
Photonic Professional 

#### Nanoscribe GmbH

Phone +49 721 981 980 0

E-Mail lifesciences@nanoscribe.com

Web www.nanoscribe.com